



NEW Adult Trainer Range

Following the introduction of our stitchdown range of paediatric footwear, Smart Steps, we have been inundated with requests to produce adult styles using the same construction, and therefore we are pleased to introduce these new trainer styles into our Chaneco Adult range. All our usual Chaneco features and benefits apply:

- Available in sizes 1 - 14 including ½ sizes
- Slim, standard and wide fitting options
- Side Stitched or Wrap-Up Trainer soles
- Supplied with 1 x 3mm poron inlay and a contoured polyurethane footbed
- Full leather linings and padded collars
- All modular and sole adaptations available

LADIES

Chelsea

DEFAULT OPTIONS:

Last - Flat

Sole - Side Stitched

COLOURS SHOWN:

Turquoise Bulbuck

Pink Leather

NHS SCHEDULE CODE:

F5c



Bailey

DEFAULT OPTIONS:

Last - Flat

Sole - Wrap-up Trainer

COLOURS SHOWN:

White Leather with Fuchsia Glitter Trim

Ice Bulbuck with Navy Leather Trim

NHS SCHEDULE CODE:

F5c



Eden

DEFAULT OPTIONS:

Last - Flat

Sole - Side Stitched

COLOURS SHOWN:

Navy Leather

Red Bulbuck

NHS SCHEDULE CODE:

F6b



Morgan

DEFAULT OPTIONS:

Last - Flat

Sole - Wrap-up Trainer

COLOURS SHOWN:

Purple Leather with White Leather Trim

Gooseberry Bulbuck with White Leather Trim

NHS SCHEDULE CODE:

F6b



FASTENINGS

All styles are available in either lace or Velcro and any combination of eyelets, hooks and straps are possible - please specify on your order



All Eyelets



Eyelets with top hook



2 Eyelets otherwise all hooks

Please note:
these styles feature stitching that may make them unsuitable for diabetic patients.

MENS

Brooklyn

DEFAULT OPTIONS:

Last - Flat

Sole - Side Stitched

COLOURS SHOWN:

Brown Leather

Red Bulbuck

NHS SCHEDULE CODE:

F5c



Campbell

DEFAULT OPTIONS:

Last - Flat

Sole - Wrap-up Trainer

COLOURS SHOWN:

Black Leather with White Leather Trim

Gooseberry Bulbuck with White Leather Trim

NHS SCHEDULE CODE:

F5c



Harley

DEFAULT OPTIONS:

Last - Flat

Sole - Side Stitched

COLOURS SHOWN:

Tan Leather

Black Bulbuck

NHS SCHEDULE CODE:

F6b



Jordan

DEFAULT OPTIONS:

Last - Flat

Sole - Wrap-up Trainer

COLOURS SHOWN:

White Leather with Black Leather Trim

Turquoise Bulbuck with White Leather Trim


NHS SCHEDULE CODE:

F6b



FLAT LAST

- Offers a pitchless plantar surface
- Round toe shape
- Standard 3mm Inlay Allowance
- Slim, standard and wide fitting options
- Restricts movement to reduce pain
- Forward propulsion without dorsiflexion

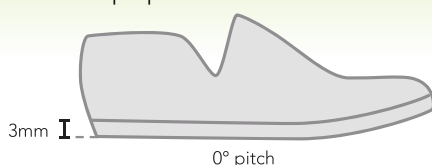

DEFAULT SOLES

All our trainer styles have a cemented construction with either a Side Stitched sole or Wrap-up Trainer Sole depending on the style in either **black or white**. Please note that the default will be supplied unless specified on the order. Both sole units have a rubber bottom which could be changed to lightweight EVA to reduce weight.



Side Stitched Sole



Wrap-up Trainer Sole

Materials & Colours

All our footwear is fully lined in pig skin leather as standard. We can also use calf leather, moleskin or sheepskin on request (additional charge).


WAXY LEATHER



BULBUCK (SUEDE)



GLITTER TRIM

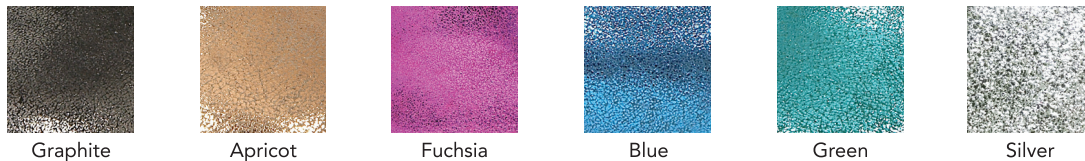

Customisation

All Chaneco footwear is made to order in Northampton by our team of experienced footwear technicians and therefore we have the ability to customise your footwear to your patient's requirements. All our other colours and trims are available for these styles on request and if you want something a bit different - just ask!



Looking After Your Footwear

Waxy Leather - Remove any surface dirt using a damp cloth. Then either use shoe polish on a soft cloth and polish to a shine or use colourless leather balm on a soft clean sponge to cover scuffs and restore to a matt finish.

Bulbuck (Suede) - Before wear, spray with suede protector. Remove any surface dirt with the suede brush or nail brush. Use white vinegar on a soft white cloth to rub any stains, allow to dry completely then brush again. Spray with suede protector.

(N.B. When footwear is new, the very strong colours in the Bulbuck (suede) material may leave some colour marking on furniture and other surfaces).

Glitter & Patent Leather - Remove any dirt using a damp cloth & use patent cleaner or petroleum jelly on a soft cloth and polish to a shine.

Lead Times

Standard and modular footwear including adaptations will be despatched within 10 working days. Bespoke footwear to a draft, scan or casts will be despatched within 15 working days. If you need footwear urgently, please phone us with your requirements and we will endeavour to meet your timescales.

Customer Services

Our Customer Services Team and your local Business Manager are all on hand to help you should you need any help regarding our range. Please call us on 01604 709999 or email info@chaneco.co.uk

Copies of our Terms & Conditions, order forms, measuring guidance and further information about Chaneco footwear can all be found on our website www.chaneco.co.uk

